




epredia

Versatile. Scalable. Flexible.

LIS Workflow Solutions and Integration Options

SlideMate AS Slide Printer
Gateway Software

Find it at fisherhealthcare.com

 **fisher**
healthcare
part of Thermo Fisher Scientific


Options to meet your workflow needs



1 SlideMate™ AS Slide Printer Label Image Workflow

PC Connected to LIS

User scans cassette or manually enters requests for slides in their LIS

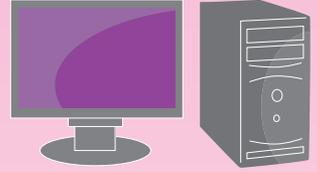

LIS system database

LIS database sends label data to SlideMate AS via print driver



1-by-1 label image data

Print driver

Sends print data to complete the job



SlideMate AS

Prints requested label data upon slide.


Discover the advantages of EpreDia SlideMate labeling systems by visiting epredia.com/slidemateAS

Have questions about **LIS integration**? Please consult with an EpreDia representative to discover which options will best support your LIS workflow.

2 SlideMate AS Slide Printer ZEBRA® Label Format Workflow


PC Connected to LIS

User scans cassette or manually enters requests for slides in their LIS

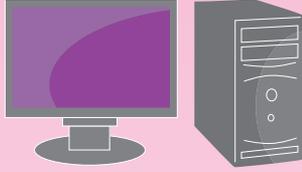

LIS system database

LIS database sends label data to SlideMate AS via print driver

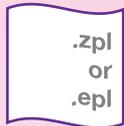

ZEBRA Label Format

Generic print driver sends data



SlideMate AS

Prints requested label data upon slide.

3 SlideMate AS Slide Printer Data String Workflow with Gateway Software

PC Connected to LIS

User scans cassette or manually enters requests for slides in their LIS



LIS system database

LIS database sends label data to SlideMate AS via print driver



Text delimited data

Shared folder



Gateway Software

Monitors the local or network share folder and moves data from the folder to the printer that will be used to create label images.



SlideMate AS

Prints requested label data upon slide.

SlideMate AS Printer Ordering Information

Systems	Order Number	LIS Compatibility
SlideMate AS Basic Printer SlideMate AS basic printer with color touchscreen, power supply, power cord and Windows® Printer Driver software to print externally generated slide label from an image file	A83900001	①
SlideMate AS Basic Printing System SlideMate AS basic printer with slide delivery system, color touchscreen, power supply, power cord and Windows® Printer Driver software to print externally generated slide label from an image file	A83900002	①
SlideMate AS On-Demand Printer SlideMate AS on-demand printer with color touch-screen, built-in scanner, power supply, power cord and advanced software enabling label creation, special panel printing and import from LIS	A83900004	① ② ③
SlideMate AS On-Demand Printing System SlideMate AS on-demand printer with slide delivery system, color touch-screen, built-in scanner, power supply, power cord and advanced software enabling label creation, special panel printing and import from LIS	A83900005	① ② ③

Accessories	Cat. Number
Thermal transfer ribbon	A83910001
Slide delivery system	A83910002
Spare slide hopper for SlideMate AS slide delivery system	A83920016
Gateway software (1 unit)	A83910005
Multi-Gateway software (up to 10 units)	A83910009
Head cleaning kit	A83910008

Please consult with an EpreDia representative to discover which options will best support your LIS workflow.

SlideMate slides and SlideMate Plus adhesion slides from EpreDia

Experience crisp, clean printing with SlideMate™ microscope slides from EpreDia, designed for the SlideMate AS slide printer*. SlideMate slides can help reduce identification issues common to handwritten labels when paired with the SlideMate AS slide printer.

*The heat setting for SlideMate AS printers when used in conjunction with SlideMate slides may need to be lowered. Each SlideMate AS printer should be optimized individually.

① **SlideMate AS Label Image Workflow**
 Select one SlideMate AS system. Windows Printer Driver included.

② **SlideMate AS ZEBRA Label (ZPL) Data Format**
 Select one SlideMate AS on-demand system.

③ **SlideMate AS Data String Workflow, Gateway Software**
 Please consult with an EpreDia representative to discover which accessory software will best support your workflow needs.



SlideMate Microscope Slide

- White tab, 25 x 75 x 1 mm, ground edges, 90° corners
 - 72 slides/pack, 20 packs/case
- US/CAN: TT-40411-8217-001
 Other: TT-40418217-W

SlideMate Plus Adhesion Microscope Slide

- White tab, 25 x 75 x 1 mm, ground edges, 90° corners
 - 72 slides/pack, 20 packs/case
- US/CAN: TT-40411-PS-8218-1
 Other: TT-40418218-PS-W

Distributed by Fisher Healthcare
 Contact us today:

In the United States

Order online: fisherhealthcare.com
 Fax an order: 1-800-290-0290
 Call customer service: 1-800-640-0640

